

5 - POINT TECHNIQUE & KEY FEATURES



5 - POINT TECHNIQUE:

1. highest point of the zygomatic bone
2. point of intersection of midpupillary line & nasal wing line
3. point in the ear tragus
4. point of the mental region (chin)
5. point in the region of the lower jaw (nerve exit point)

INTRODUCTION

Kamilane Elite Timeless Touch

Whitening & Anti-Wrinkle / 3ml x 5 syringes

Kamilane Elite offers a **comprehensive and advanced approach** to skin rejuvenation, designed to **restore and maintain** a radiant, youthful, and healthy complexion. Through regular treatments, patients can experience **long-term improvements** in skin health, hydration, and overall vitality.

This innovative solution is **ideal for individuals** seeking to enhance skin **hydration, elasticity, and firmness**, promoting a naturally revitalized appearance. It is particularly **beneficial for those looking to diminish** the visibility of **fine lines, wrinkles, and other signs of aging**, resulting in smoother and more youthful-looking skin. Additionally, Kamilane Elite is well-suited for **patients interested in brightening** their complexion, correcting uneven skin tone, and achieving a **more luminous and balanced** skin texture.

Collagen, Elastin, and Hyaluronic Acid (HA) are three key components essential for maintaining beautiful skin. Collagen provides **firmness**, Elastin ensures **elasticity**, and HA supports **regeneration**.


NANO-MESH CROSS-LINKING

Kamilane Elite utilizes advanced Nano Mesh Cross-linking (NMC) technology to safely replenish the skin's HA levels. Thanks to the product's molecular stability, its benefits are long-lasting and become increasingly evident over time.

WHITENING & ANTI-WRINKLE

Kamilane Elite is available in three advanced formulations: Anti-Wrinkle, Whitening, and a combined Whitening & Anti-Wrinkle solution. It is designed to reduce fine lines, wrinkles, and signs of aging, promoting a youthful look.

ACTIVE INGREDIENTS

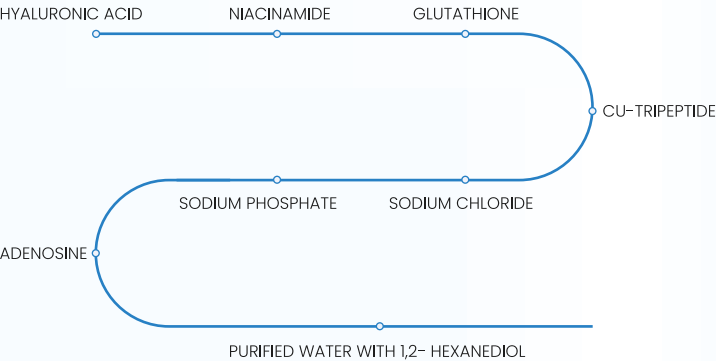


Kamilane

Timeless Touch
ELITE

WHITENING & ANTI-WRINKLE
3ml x 5 syringes

PREMIUM QUALITY

Kamilane Elite is a premium injectable skin booster engineered to address multiple visible signs of aging, such as wrinkles, fine lines, and skin dehydration. This highly effective treatment, formulated with a potent blend of advanced ingredients, enhances hydration, boosts elasticity, and promotes a brighter, more even complexion. Known for its exceptional rejuvenating benefits, Kamilane Elite is a popular choice in South Korea for achieving radiant and youthful skin.

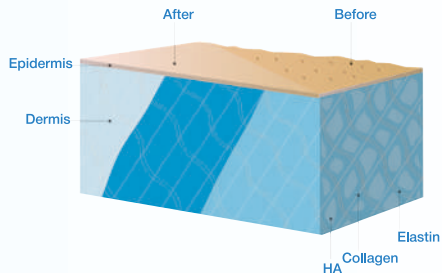

KEY BENEFITS



Improved skin elasticity:
Boosts collagen production, enhancing the skin's resilience and firmness.

Comprehensive rejuvenation:
Improves skin texture, tone, and firmness with continued use, for a visibly refreshed appearance.

Skin brightening & Tone correction:
Hydroquinone and glutathione lighten pigmentation and even skin tone, revealing a radiant complexion.

Collagen stimulation:
Cur-Tripeptide enhances collagen and elastin production, restoring skin firmness and elasticity.

Long-lasting anti-aging effect:
Reduces fine lines and wrinkles, improving overall skin smoothness and youthful texture.

Deep hydration:
Provides lasting moisture by attracting and retaining water in the skin, combating dryness and maintaining a plump, smooth appearance.

Patient comfort
Constructed with high-quality materials for a safe, smooth, and reliable user experience.

Safe application
Engineered and rigorously tested by experts to ensure precision and safety.

Innovative design
The conical plunger stopper enhances uniformity and smooth application for optimal performance.

Comfortable handling
Constructed with high-quality materials for a safe, smooth, and reliable user experience.

Enhanced visual orientation
Ultra-thin wall needles minimize bruising and reduce side effects, ensuring patient comfort.

APPLICATION PROTOCOL

5 - Point Technique Protocol

- 1

5 - POINT TECHNIQUE INJECTION POINTS

Highest point of the zygomatic bone
Intersection of midpupillary line & nasal wing line
Ear tragus
Mental region (chin)
Lower jaw (nerve exit point)
- 2

FACE INJECTION PROTOCOL

Volume: 2ml
Technique: 5 injection points per side (1ml per side)
Injection Volume: 0.2ml per point
Depth: Subdermal layer
- 3

NECK INJECTION PROTOCOL

Volume: 2ml
Technique: 10 injection points across the neck
Injection Volume: 0.2ml per point
Depth: Preplatysmal fat
Coverage: Even distribution across the neck for optimal rejuvenation
- 4

COURSE OF PROCEDURES AND RESULTS

Total Treatments: 3 procedures (2.2ml each time)
Frequency: 1 session every 14 days
Results: The patient will notice lightening and increased turgor within 24 hours. The final result can be assessed strictly after 28 days.

Microinjections Technique Protocol

- 1

INJECTION PROTOCOL

Administer via intradermal microinjections, ensuring precise epidermal-level delivery.
Utilize a multi-needle injection method for uniform diffusion and enhanced efficacy.
Space injections 1 cm apart, delivering 0.2 ml per point, increasing to 0.4 ml for deeper wrinkles.
- 2

RECOMMENDED DEVICES

Microneedling (1–4mm needles)
Hydra Stamp / Pen / Roller
EZ Injector
- 3

COURSE OF PROCEDURES AND RESULTS

Initial Session: Establishes baseline skin revitalization.
Follow-Up: Conducted 4–6 weeks post-treatment for sustained results.
Maintenance: Adjusted based on individual response and clinical assessment.
Results: Glowing skin appears within 3–5 days post-injection. Effects last approximately 4 weeks, with continued improvement through regular treatments.

Post-Treatment Care & General Precautions

Advise patients to avoid excessive heat, direct sun exposure, and harsh skincare for 48–72 hours.

For Professional Use Only: To be administered by licensed medical professionals trained in injectable procedures.



KS BIO
CONTACT

Address: #301, 3rd floor, Kim's Tower 2, 637-6 Myeonmok-dong, Jungnang-gu, Seoul, Republic of Korea
E-mail: info@ksbiokorea.com
Website: www.ksbiokorea.com